

Connecting PostureScribe to ClinicMind

The one-time practice setup, start to finish.

i Already using PostureScreen or PostureRay? PostureScribe uses the same authentication mechanism those products already use to post X-rays and images into ClinicMind. If the practice is connected, Scribe documents just work once the feature is enabled in the app.



1 Section 1 • One-Time Practice Setup In ClinicMind (web) · once per practice, not per practitioner or patient

PART A • GET THE API KEY

Generate or locate the API key

In the **left-hand navigation's search bar**, type **"patient file upload"** and open **Patient File Upload Service**.

Under *Select the Vendor*, choose **postureScreen**, click **Activate Service**, then copy the **Key**.

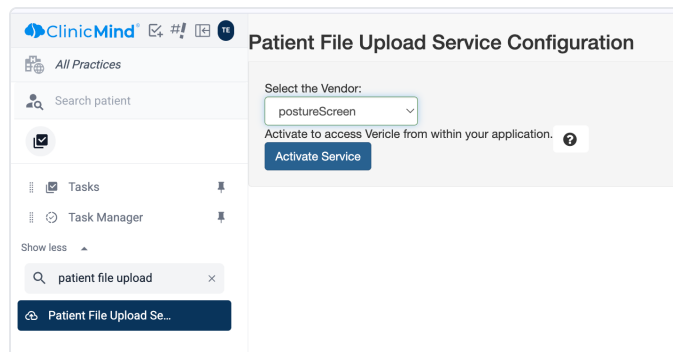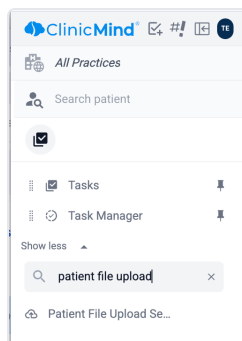


⚠ PostureScreen and PostureRay each use their own key — don't reuse one key across both.

PART B • GET THE PHYSICIAN ID

Find the Physician ID

In the **left-hand navigation's search bar**, type **"settings"** and open **Settings**.

On the **Practitioner** tab, copy the value in the **Pract. ID** column. Only one Physician ID is needed per practice.

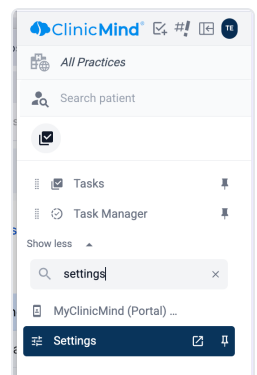

PROV	Practitioner	SvcFac	PinMap	Log	CP CPTS	Settings	Invoice	POS	P
Pract. ID	Person ID	Last Name	First Name	M.I.	Schedule Name	NPI			
JOEGREEN3		GREEN	JOE		SHOCKWAVE THERAPY				
JOEGREEN2		GREEN	JOW		DECOMPRESSION				
JOEGREEN1		GREEN	JOE	A	DR. GREEN	4795178320			

⚠ Label it exactly "Physician ID" — lowercase "d". "Provider ID" or other capitalization causes confusion.

2 Section 2 • Connect the App On the practitioner's iPad

Open the PostureScreen app on the iPad and log in to the existing PostureScreen account.

✓ Existing integration: the key is likely already saved — no copy/paste needed.

📁 Brand-new practice: paste the API key from Part A into the app.

REMEMBER

Setup is **once per practice** — not per practitioner or patient. One Physician ID. Separate API keys per product. After that, PostureScribe documents post into the patient file automatically.